


MTC Bio™ DISPOSABLE SPREADERS

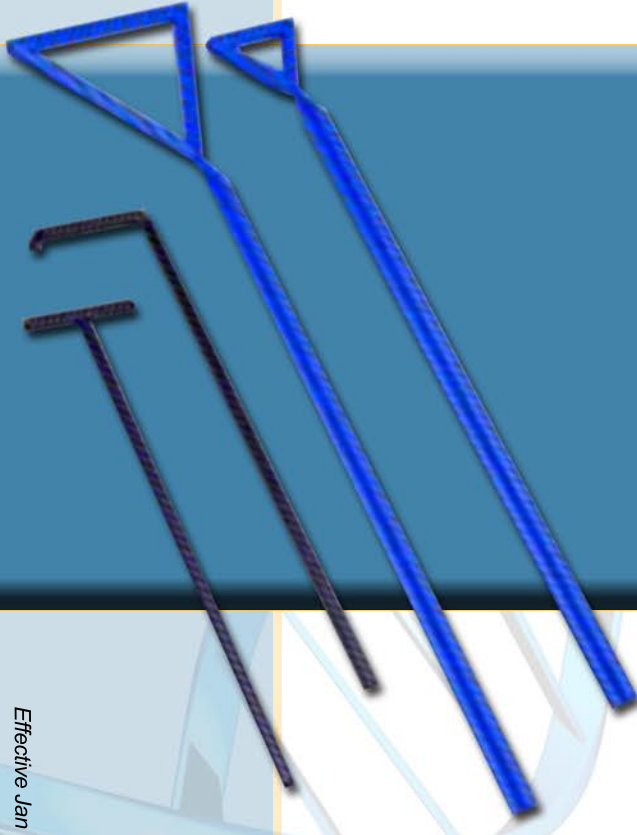

- **Spread cultures & liquids evenly**
- **Smooth-edges - minimize “gouging”**
- **3 Configurations: “L”, “T”, or Triangle**
- **Ideal for bacterial CFU counts**
- **Presterilized - no flaming required**

The MTC Bio line of bacteriological spreaders has been expanded to include four different shapes and sizes. In addition to the original L-shape, a similar T-shape version and two triangle shapes (large and small) are now available. All versions are supplied in sterile packaging and have smooth, rounded edges to prevent unwanted cutting, scraping or gouging of agar during inoculation.

L-shape The original spreaders are ideal for spreading bacteriological cultures evenly. The shape allows the user to rotate and apply easily throughout an entire round plate or dish. They are made of ABS and have a slight upward curve at the edge to prevent dragging or scraping.

T-shape This shape ensures that consistently even pressure can be applied across the entire span of the 34mm spreading surface. Designed for spreading and dispersing liquids across the surface of agar culture plates, they are made of ABS and commonly used for CFU plate counts on water, milk, urine, etc.

Triangle Shape This new design is also used for CFU count applications, but has a grooved handle that improves grip and orientation. Made of polypropylene, they are available with a 30mm span or an extra large 60mm span.



2600A Main Street Extension
Sayreville, NJ 08872
Ph: 908-279-0280 www.mtcbio.com

Effective Jan 1, 2017. Subject to change without notice.

Ordering Information

Cat. No.	Description	Dimensions	Sterile Unit	Units/Pack	Qty/Case	Price
M5421	Spreader, L-shape	38 x 156mm	10 per peel pack	50 peels	500	\$77.00
M5420	Spreader, L-shape	38 x 156mm	1 per peel	500 peels	500	\$89.00
M5431	Spreader, T-shape	34 x 140mm	10 per peel pack	50 peels	500	\$85.00
M5430	Spreader, T-shape	34 x 140mm	1 per peel	500 peels	500	\$95.00
M5441	Spreader, Triangle	60 x 235mm	10 per peel pack	10 peels	100	\$59.00
M5440	Spreader, Triangle	60 x 235mm	1 per peel	500 peels	500	\$320.00
M5451	Spreader, Triangle	30 x 207mm	10 per peel pack	10 peels	100	\$59.00