

SureStrain™ Premium Cell Strainers

- Advanced design for faster flow through
- Interlocking and stackable design
- Designed for use with 50mL tubes
- RNase, DNase, Pyrogen Free
- Individual sterile blister packaging




The SureStrain™ Premium Cell Strainer system is a fast and easy-to-use solution for detaching, isolating and straining cells from clusters and tissues. A polypropylene frame supports a nylon mesh that is available in three color-coded porosities: 40µm (blue), 70µm (white), and 100µm (yellow). The advanced frame/housing design also contributes to superior performance and reduced processing time, by eliminating the “air-lock” effect. Improved aeration allows air to escape while filtrate is entering the tube, preventing messy spillovers.



Even more features to see in our video demonstration!

- Reducing adapter allows coupling with smaller tubes
- Elevated feet enable straining into flat dishes
- Molded rim for suspension in 6-well plates

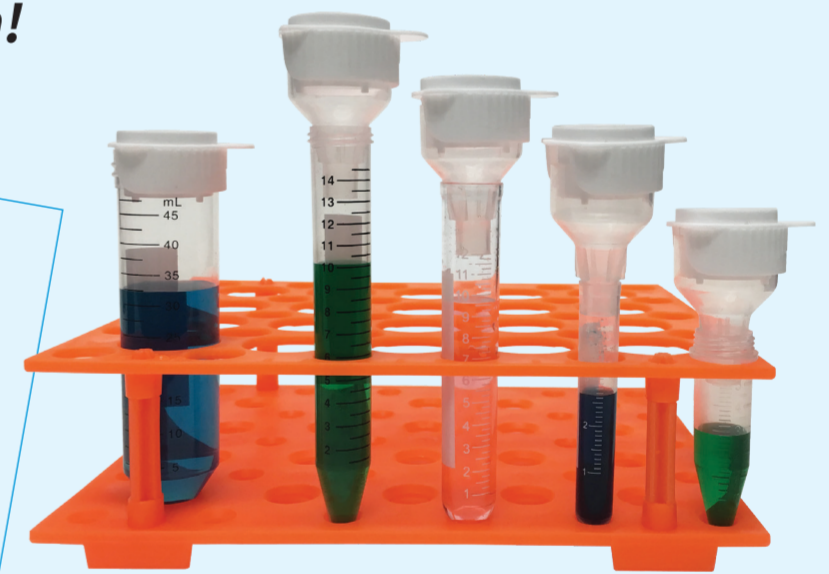

CAT. NO.	DESCRIPTION	QTY	LIST PRICE
C4040	40µm Cell Strainer, Blue, w/ 1 Reducing Adapter	50/pack	\$109.00
C4070	70µm Cell Strainer, White, w/ 1 Reducing Adapter	50/pack	\$109.00
C4100	100µm Cell Strainer, Yellow, w/ 1 Reducing Adapter	50/pack	\$109.00
C4010	PolyPestle™ Double-Sided Pestle, Individually Wrapped - for use with cell strainers or microcentrifuge tubes	50/pack	\$51.00
C4015	Reducing Adapter, Multi-Width Funnel for Smaller Tube Sizes - for use with 15mL, 17×100mm, 12×75mm, and 5mL tubes	25/pack	\$47.88